

Softlift™ Laparotomy vacuum manipulator

The Softlift™ is a single use, non-invasive, vacuum activated all-in-one tissue manipulation device.


The Softlift™ can be used in a laparotomy, mini-laparotomy and transvaginal procedures by grasping, retracting and manipulating internal organs/tissues atraumatically.

A small vacuum cup made of soft, deformable material allows for easy cup insertion into narrow cavities. The cups flat profile and flexible stem enables generous internal visualisation therefore quicker completion of the procedure.

The cup has vertical ridges branching from the centre with attachment to the tubing at a 45° angle which provides greater traction quality with minimised probability of slipping along the tissue.

The soft material allows for the penetration of an aspiration needle through the cup without loss of vacuum.

The PalmPump™ provides a safe and effective source of vacuum right in the palm of your hands. The operator has full control and the unique traction device allows flexibility in manipulation.

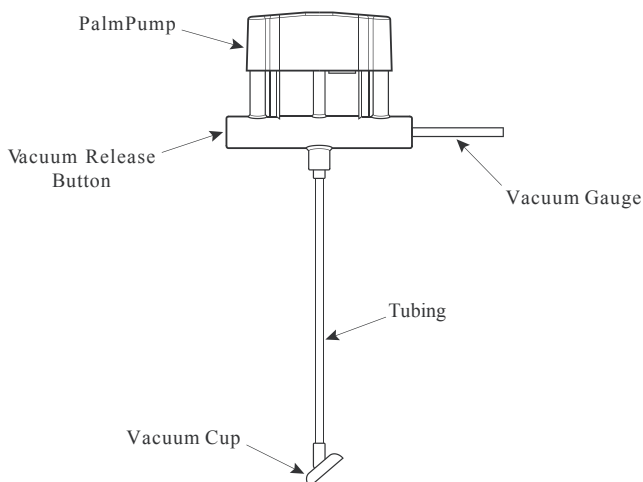


Ordering Information

Part Number	LVM-1000
Quantity	5 units per case
Delivery	Usually ex-stock
Description	Softlift™

Softlift™ can be used in a variety of procedures, including:

Tubal ligations
Hysterectomies
Oophorectomies
Myomectomies
Ectopic pregnancy
Cystectomies
Cholecystectomies
Tuboplasties
Dermoid cysts
Endometriomas
Transvaginal hysterectomies
Bowel

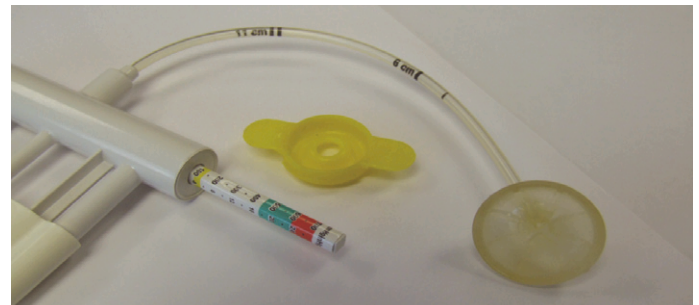

www.clinicalinnovationseurope.com

Distributed